


EMTEK Microbial Air Sampler Comparison Data

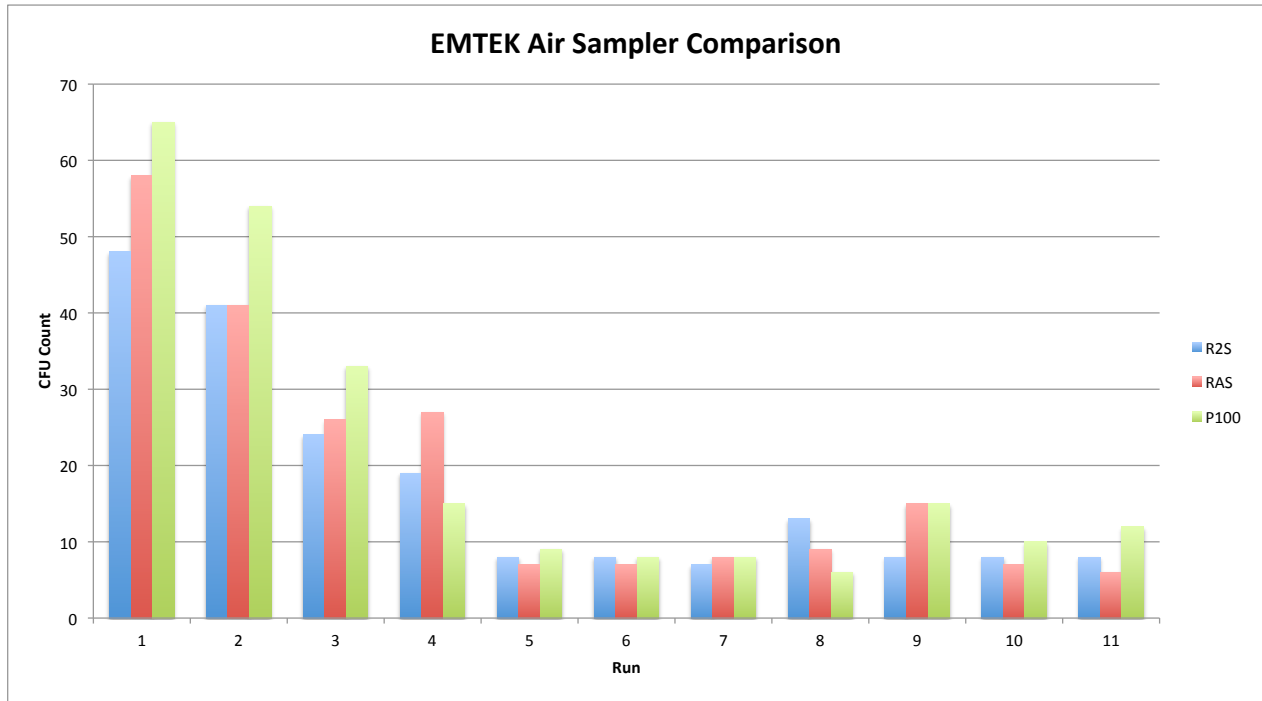

TEST DATE:		9-Aug-13	Final Read - 13AUG2013 (13:15)			
START TIME	VOLUME (Liters)	RUN	CFU/Run			
			R2S	RAS	P100	
10:25	1000	1	48	58	65	
11:08	1000	2	41	41	54	
11:50	500	3	24	26	33	
12:11	500	4	19	27	15	
12:37	250	5	8	7	9	
13:23	250	6	8	7	8	
13:35	250	7	7	8	8	
13:49	250	8	13	9	6	
14:05	250	9	8	15	15	
14:17	250	10	8	7	10	
14:29	250	11	8	6	12	
			R2S	RAS	P100	
Total Volume	4750	Total CFU	192	211	235	
			CFU/L	0.040	0.044	0.049
			CFU/1000L	40.42	44.42	49.47

NOTE: Flow rate of 28.3 LPM used with all devices.